KC-SMART Vignette

Jorma de Ronde, Christiaan Klijn April 22, 2010

1 Introduction

In this Vignette we demonstrate the usage of the KCsmart R package. This R implementation is based on the original implementation coded in Matlab and first published in January 2008 [1].

KC-SMART is a method for finding recurrent gains and losses from a set of tumor samples measured on an array Comparative Genome Hybridization (aCGH) platform. Instead of single tumor aberration calling, and subsequent minimal common region definition, KC-SMART takes an approach based on the continuous raw log2 values

Briefly: we apply Gaussian locally weighted regression to the summed totals positive and negative log2 ratios separately. This result is corrected for uneven probe spaceing as present on the given array. Peaks in the resulting normalized KC score are tested against a randomly permuted background. Regions of significant recurrent gains and losses are determined using the resulting, multiple testing corrected, significance threshold.

2 Basic Functionality

After installing the KCsmart package via Bioconductor, load the KCsmart code into R.

> library(KCsmart)

Load some sample aCGH data. We supply an example dataset of an old 3000 probe BAC clone experiment. KC-SMART has no problem running on modern 300k+ probe datasets, but this dataset is supplied to keep downloading times at a minimum. This dataset can be a simple dataframe, DNAcopy object or a CGHbase object, and not contain any missing values. If the data contains missing values, we advise the user to either impute the missing values using the mean of the two neighboring probes or to discard probes containing missing values.

> data(hsSampleData)

This data is the simplest form available to use in KCSMART. As can be observed using:

> str(hsSampleData)

```
'data.frame':
                    3268 obs. of 22 variables:
                  "1" "1" "1" "1" ...
$ chrom
         : chr
$ maploc
           : num
                  200000 200050 3284807 4542451 7103443 ...
$ sample 1 : num
                  1.326 0.0368 0.4335 1.8879 1.2002 ...
$ sample 2 : num
                  1.1155 -0.0478 0.641 0.456 0.831 ...
$ sample 3 : num
                 0.4204 0.6486 -0.3638 -0.2573 0.0953 ...
$ sample 4 : num
                  0.06313 -0.00108 0.13552 0.57361 -0.09144 ...
$ sample 5 : num
                  0.816 0.843 0.893 1.275 0.686 ...
$ sample 6 : num
                  0.471 1.262 0.4 0.987 0.725 ...
  sample 7 : num
                  1.296 0.471 0.614 1.007 1.041 ...
                  0.835 0.522 0.643 0.572 0.953 ...
  sample 8 : num
                  1.98 -0.529 0.808 -0.165 1.668 ...
  sample 9 : num
$ sample 10: num
                  0.668 0.313 1.083 0.7 1.782 ...
  sample 11: num
                  0.338 -0.118 0.205 1.402 0.299 ...
$ sample 12: num
                  0.688 -0.187 0.593 0.941 0.625 ...
                  0.4127 0.5349 0.9524 0.9948 0.0954 ...
$ sample 13: num
                  0.938 0.254 0.599 0.615 0.466 ...
$ sample 14: num
  sample 15: num
                  0.0481 -0.141 0.9038 0.5821 0.2582 ...
$ sample 16: num
                  0.4 -0.306 0.198 0.652 0.589 ...
                  1.178 0.72 0.731 1.192 0.287 ...
$ sample 17: num
$ sample 18: num
                  0.2878 0.5926 1.3355 0.6192 -0.0581 ...
$ sample 19: num
                  0.738 0.489 0.739 1.075 -0.109 ...
$ sample 20: num
                 1.552 0.439 0.264 0.841 0.218 ...
```

The user can easily model their data in this format, basically only requiring a column named \$chrom containing the chromosome on which the probe is located and a column named \$maploc containing the chromosomal location of the probe. Alternatively, both DNAcopy and CGHbase objects can be used as direct input for the calcSpm function discussed below.

KC-SMART uses mirroring of the probes at the ends of chromosomes and near centromeres to prevent signal decay by the convolution. The locations of the centromeres and the lengths of the chromosomes are therefore necessary information. Hard-coded information about the human and mouse genome is supplied, but the user can easily provide the coordinates themselves in R. The 'mirrorLocs' object is a list with vectors containing the start, centromere (optional) and end of each chromosome as the list elements. Additionally it should contain an attribute 'chromNames' listing the chromosome names of each respective list element.

To load the presupplied information use:


> data(hsMirrorLocs)


For an analysis in mouse data(mmMirrorLocs) is used. Using the calcSpm command the convolution is performed at the given parameters and stored in a samplepoint matrix. Here we perform the convolution at two different kernel-widths (first the default sigma = 1Mb, and a second 4Mb sigma). We apply the default parameters to run the convolution.

```
> spm1mb <- calcSpm(hsSampleData, hsMirrorLocs)
> spm4mb <- calcSpm(hsSampleData, hsMirrorLocs, sigma = 4e+06)</pre>
```

Now we plot the chromosome-wide view of the convolution for both the gains and the losses for the 1Mb kernelwidth analysis. Since the samplepoint matrix we just created is an S4 object we can just apply the plot command to it to visualize the results.

> plot(spm1mb)

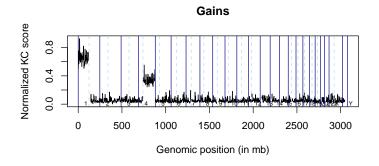


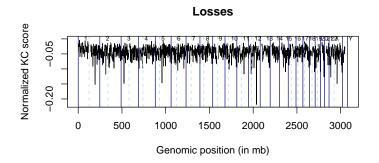

It is possible to select only particular chromosomes in the plot command, and to only display the gains or only the losses. We will now visualize the KCsmart output for gains of chromosomes 1, 12 and X.

> plot(spm1mb, chromosomes = c(1, 12, "X"), type = "g")

Gains

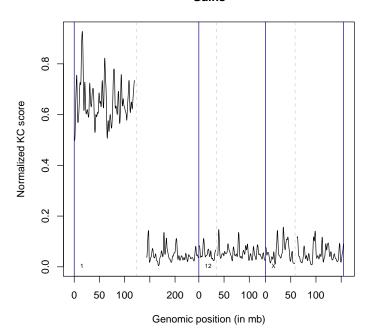

As can be seen, in this particular (synthetic) dataset chromosome 1 shows a large recurrent gain compared to other chromosomes. How significant is this gain? We can determine a significance level based on permutation to determine this. In this example we will permute the data 10 times and determine a threshold at p < 0.05. We recommend 1000 permutations if you want to do a definitive analysis of a dataset.

```
> sigLevel1mb <- findSigLevelTrad(hsSampleData, spm1mb, n = 10, + p = 0.05)

[1] "Calculating alpha = 0.05 significance cut-off"

[1] "Found 553 pos peaks and 567 neg peaks in observed sample point matrix"


[1] "Starting permutations .."

At iteration 1 of 10
At iteration 2 of 10
At iteration 3 of 10
At iteration 4 of 10
At iteration 5 of 10
At iteration 6 of 10
At iteration 7 of 10
At iteration 8 of 10
At iteration 9 of 10
At iteration 9 of 10
At iteration 10 of 10[1] "\n"
```

We can now plot the significance level we just calculated in the genomewide view of the KCsmart result using the sigLevels= parameter of the plot function. This time, we'll use the type=1 parameter. This will plot the gains and losses in a single frame.

> plot(spm1mb, sigLevels = sigLevel1mb, type = 1)

Gains and losses

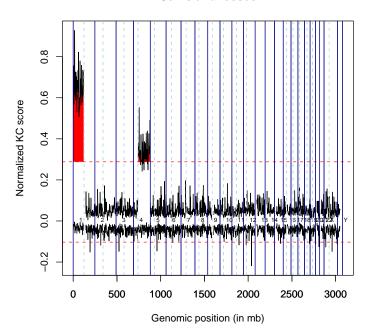

Let's take a closer look at the gains of chromosome 1, 12 and X again, now with the added significance level.

> plot(spm1mb, chromosomes = c(1, 12, "X"), type = "g", sigLevels = sigLevel1mb)

Gains

We've also implemented a less stringently corrected version of the permutation algorithm that was shown above. This approach doesn't calculate the Bonferroni corrected significance level, but the False Discovery Rate (FDR). The command to find the FDR for this particular dataset would be FDR1mb <- find-SigLevelFdr(hsSampleData, spm1mb, n=10, fdrTarget = 0.05). In general the FDR is less conservative, and caution must be taken when using this as a threshold.kile

100 0

50 100

50

Genomic position (in mb)

Now we have our significance level we would like to get usable information back about which regions are found to be significantly abberrant and which probes of the original dataset are contained in these regions. Using the get—SigSegments function we can input a samplepoint matrix and our calculated significance level to get back the relevant information.

> sigRegions1mb <- getSigSegments(spm1mb, sigLevel1mb)</pre>

Now let's have a look at the information contained in sigRegions1mb, it is again an S4 object and will display its contents if you just type its name.

> sigRegions1mb

Significantly gained and lost segments Cutoff value was 0.288005792169466 for gains and -0.103074323665016 for losses A total of 10 gains and 15 losses were detected

Use 'write.table' to write the file to disk

50 100

0

You could write the infomation contained in sigRegions1mb to a file using write.table(write.table(sigRegions1mb, file='sig.txt'). You can also

access the information in R using subsetting. For example, if you wanted to find the probe-names of the probes contained in the first significantly gained region you can get them like this:

> sigRegions1mb@gains[[1]]\$probes

```
2
                       4
                            5
  [1]
                  3
                                     7
                                          8
                                               9
                                                  10
                                                       11
                                                            12
                                                                13
                                                                     14
                                                                          15
                                                                               16
                                                                                   17
                                                                                        18
 [19]
       19
            20
                 21
                      22
                          23
                               24
                                    25
                                         26
                                             27
                                                  28
                                                       29
                                                            30
                                                                31
                                                                     32
                                                                          33
                                                                               34
                                                                                   35
                                                                                        36
 [37]
       37
            38
                 39
                      40
                          41
                               42
                                    43
                                         44
                                             45
                                                  46
                                                       47
                                                            48
                                                                49
                                                                     50
                                                                          51
                                                                               52
                                                                                   53
                                                                                        54
                               60
                                                                                   71
                                                                                        72
 [55]
       55
            56
                 57
                      58
                          59
                                    61
                                         62
                                             63
                                                  64
                                                       65
                                                            66
                                                                67
                                                                     68
                                                                          69
                                                                               70
                          77
 [73]
       73
            74
                 75
                      76
                               78
                                    79
                                         80
                                             81
                                                  82
                                                       83
                                                            84
                                                                85
                                                                     86
                                                                          87
                                                                               88
                                                                                   89
                                                                                        90
 [91]
       91
            92
                 93
                      94
                          95
                               96
                                    97
                                         98
                                             99
                                                 100 101 102
                                                               103 104
                                                                        105
                                                                             106
                                                                                  107
                                                                                      108
[109] 109 110 111 112 113 114 115 116 117 118 119 120 121 122 123 124 125 126
[127] 127 128 129
```

Using str(sig.regions.1mb) you can find out more about the way the significant regions are stored in the sigRegions object.

3 Multi-scale analysis

One of the powerful features of KC-SMART is the ability to perform multi-scale analysis. The algorithm can be run using different kernel widths, where small kernel widths will allow you to detect small regions of aberration and vice versa using large kernel widths large regions of aberration can be detected. Remember our convolution step? We calculated the convolution at two sigma's: 1Mb and 4Mb. We will now calculate a significance level for the 4Mb samplepoint matrix and plot the results of both kernelwidths in a scale space figure.

```
> sigLevel4mb <- findSigLevelTrad(hsSampleData, spm4mb, n = 10)
```

[1] "Calculating alpha = 0.05 significance cut-off"

```
[1] "Found 171 pos peaks and 181 neg peaks in observed sample point matrix"
[1] "Starting permutations .."

At iteration 1 of 10


At iteration 2 of 10
```

```
At iteration 1 of 10
At iteration 2 of 10
At iteration 3 of 10
At iteration 4 of 10
At iteration 5 of 10
At iteration 6 of 10
At iteration 7 of 10
At iteration 8 of 10
At iteration 9 of 10
At iteration 10 of 10[1] "\n"
```

Now you can plot the scalespace of the two different kernel widths using plotScaleSpace. A scale space can show you small and large significant regions in one plot, as well as give you insight within a single significant region. Darker red means more significant. As an option you can plot either the gains or the losses or both. When plotting both, two x11 devices will be opened. Using the 'chromosomes' parameter the chromosomes to be plotted can be set.

```
> plotScaleSpace(list(spm1mb, spm4mb), list(sigLevel1mb, sigLevel4mb),
+ type = "g")
```



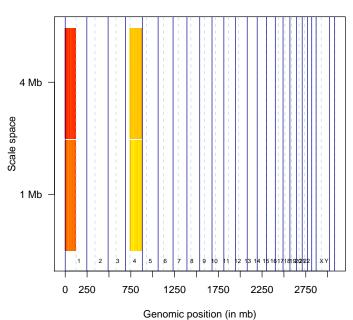

4 Comparative KC-SMART

Using KC-SMART you can also find regions of significant different copy number change between two specific groups of samples. The resulting regions of significant difference are not required to be significantly recurrent in the entire tumor set. We employ kernel smoothing on a single tumor basis, followed by either a permutation based significance analysis using the signal to noise ratio of the sample point matrices or the use of the siggenes package to determine significance.

The main function of the comparative version of KC-SMART is calcSpm-Collection, this function requires the data to be in the same form as needed by calcSpm, but it also requires a class vector. This is a vector of ones and zeros, with as many elements as there are samples in the dataset. The ones denote one class, and the zeros the other.

Our sample data contains an amplification on chromosome 4 that is specific to a subset of the samples, namely the last 10 samples. We will now calculate the significant difference between the first 10 samples and the last 10 samples in the example dataset using comparative KC-SMART

First we'll want to calculate a sample point matrix collection, the starting point for a comparative analysis. The class vector we use is just 10 zeros followed by 10 ones. Note that this procedure can take a long time if you have a large dataset! Make sure you do it only once. If you have multiple subgroups within

a dataset you can define alternative class vectors later in the analysis. You can perform this analysis again on different scales, here we use a sigma of 1 Mb.

```
> spmc1mb <- calcSpmCollection(hsSampleData, hsMirrorLocs, cl = c(rep(0,</pre>
      10), rep(1, 10)), sigma = 1e+06)
[1] "Mirror locations looking fine"
Processing sample 1 / 20
Processing sample 2 / 20
Processing sample 3 / 20
Processing sample 4 / 20
Processing sample 5 / 20
Processing sample 6 / 20
Processing sample 7 / 20
Processing sample 8 / 20
Processing sample 9 / 20
Processing sample 10 / 20
Processing sample 11 / 20
Processing sample 12 / 20
Processing sample 13 / 20
Processing sample 14 / 20
Processing sample 15 / 20
Processing sample 16 / 20
Processing sample 17 / 20
Processing sample 18 / 20
Processing sample 19 / 20
Processing sample 20 / 20
```

We get a few warnings that some sample points are listed as NA. These sample points are the mirrored sample points that fall outside the chromosome boundaries, so they will not be used in the analysis. The next step is to calculate the significance of the diffences between the subgroups. We use the compare-SpmCollection command for this. Using a parameter called method we define the method to determine significance. The option "perm" will use class label permutation to define an empirical null-distribution against which the actual differences are measured. The "siggenes" option will utilise the siggenes package to determine the significance, this option is both faster and uses less memory.

```
> spmc1mb.sig <- compareSpmCollection(spmc1mb, nperms = 3, method = c("siggenes"))
> spmc1mb.sig

Comparison of 20 samples (10 vs. 10)
Using siggenes to find significant regions
```

Now we can use the getSigRegionsCompKC command to extract the significant regions from the original data spmc1mb using the significance information contained in spmc1mb.sig.


```
> spmc1mb.sig.regions <- getSigRegionsCompKC(spmc1mb.sig)
> spmc1mb.sig.regions
```

Significantly different regions using siggenes delta value was 2.520453 , FDR = 0.01

	startrow	${\tt endrow}$	chromosom	e st	tartposition	${\tt endposition}$
1	14843	15758	•	4	52450001	98200001
2	15775	16280	•	4	99050001	124300001
3	16302	17182		4	125400001	169400001
4	17207	17605		4	170650001	190550001
5	21654	21655		6	20900001	20950001
Us	e 'write	.table'	to write	the	file to dis	K

The information contained in <code>spmc1mb.sig.regions</code> variable can be saved to disk using the <code>write.table</code> command if further processing in other programs is required. An overview of the comparative analysis can also be plotted using the <code>plot</code> function.

> plot(spmc1mb.sig, sigRegions = spmc1mb.sig.regions)

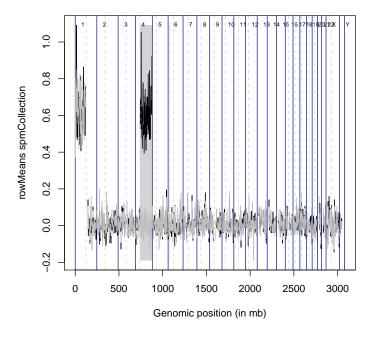

5 References

1. Klijn et.~al. Identification of cancer genes using a statistical framework for multiexperiment analysis of nondiscretized array CGH data. Nucleic~Acids~Research.~36~No.~2,~e13.